



## EQUINE NASOGASTRIC SIMULATOR

The VSI Equine Nasogastric Simulator is an integral part of clinical studies for veterinary practice. Our equine simulator is realistic in presentation and size and provides students the opportunity to use anatomical features to practice clinical procedures in a controlled setting. Through simulation students will gain confidence and expertise prior to moving to live animals.

Additionally, it is compatible with the Equine Dystocia Simulator.

## DAY ONE SKILLS & PRODUCT BENEFITS

### Nasogastric Intubation

- Replica nasal passageway that simulates anatomical structure while students practice passing the intubation tube via the esophagus.
- Well-rendered internal anatomy will allow the student to confirm tube placement in the esophagus via a visual reference site on the left side of the neck. Students can then continue to pass the tube into the stomach chamber.
- Smell and/or suction test can be performed to indicate proper placement of the tube in the stomach chamber.
- Lavage technique and administration of fluids can be simulated.

### Venipuncture and Intramuscular Injection

- Jugular venipuncture with palpable veins
- Easily replaceable latex vein tubing that can endure multiple punctures for simulation purposes
- Mechanical vein pressurizing system allows for positive pressure for blood draw (20-22 gauge needle recommended for simulation).
- Intramuscular site that can endure multiple punctures and accept fluid injections
- Soft silicone head with flexible ears and well rendered anatomy enables students to locate facial injection sites. Head also accepts bit & bridle.
- Rolling stand, adjustable for height, rostral tilt and 90 degree rotation

Equine Nasogastric Simulator with Rolling Stand 1025594

Equine Nasogastric Simulator without Rolling Stand 1025593




### TECHNICAL SPECIFICATIONS

- VSI metal rolling stand
- Durable silicone with metal armature
- Removable head for ease of transportation, and replacement of components
- Removable/washable hide covering. Can be replaced when worn
- Comes with 5 sets of replacement veins

### REPLACEABLE PARTS

- Veins 1025577
- Intramuscular Pads 1025600
- Mane 1025571
- Hide Covering 1025569

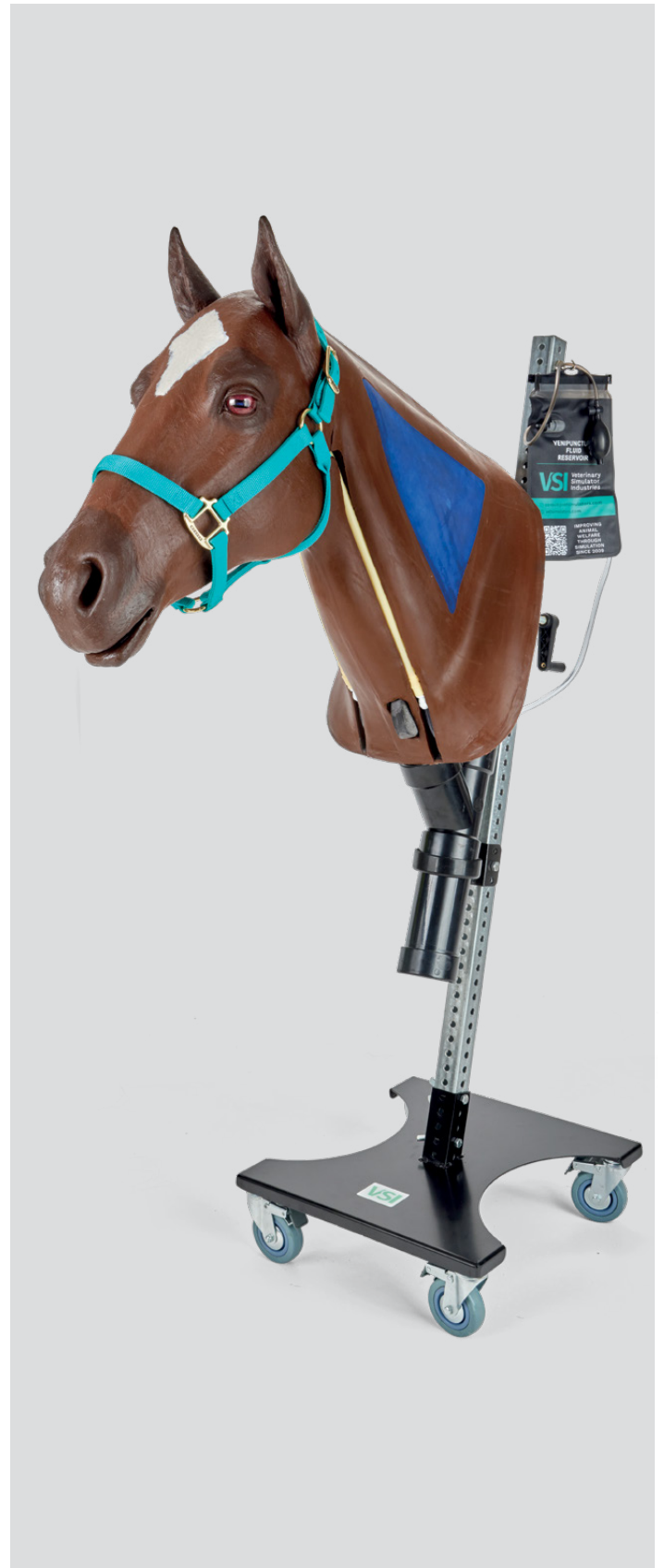

03/2025

IF YOU HAVE ANY QUESTIONS OR NEED SUPPORT, PLEASE SCAN HERE



3B Scientific Group

Headquarters · Ludwig-Erhard-Strasse 20 · 20459 Hamburg · Germany · [3bscientific.com](http://3bscientific.com)